

PHILIPS

Therapy Mask System

3100 NC
3100 SP

Nasal Cradle
Silicone Pillows

**Everything clicks,
right from the start**



The mask system that puts comfort in front

From its small size to its innovative design, Philips Therapy Mask System 3100 NC and 3100 SP with adjustable headgear is designed to help your patients get the right fit every night.^{1,2,3}

It is the only⁷ tube-in-front CPAP mask system that provides the ability to switch between a nasal cradle under-the-nose and silicone pillows cushion on the same frame. Easy to use, quiet and comfortable to wear, Philips Therapy Mask 3100 NC and 3100 SP are the perfect pair for a good night's sleep.^{1,2,3}



Easy to use

Innovative, polarized magnets in the cushion and frame allow effortless assembly of the cushions to the frame.⁴ When patients feel the magnets click into place, they can be confident their cushion is correctly attached. The intuitive design with orientation markers on the headgear for placement on the head can minimize patient mask errors that lead to callbacks, saving you time.



A comfortable fit

The majority of patients studied said that the 3100 NC and 3100 SP adjustable headgear allows for the right mask fit every night.^{1,2} The innovative tube-in-front CPAP mask system gives patients the ability to easily swap between a nasal cradle under-the-nose or silicone pillows cushion while keeping the same frame, so patients can choose what feels right for them.^{1,2,4} 3100 NC and 3100 SP's quiet sound, small size and light weight are designed to be less disruptive for patients.^{5,6}



A great value

Philips Therapy Mask System 3100 NC and 3100 SP offer high value by creating workflow efficiencies with its unique features that can reduce setup time as well as the need for high-touch refits. Having the choice of nasal cradle under-the-nose and silicone pillows cushions that use the same sizing gauge and the same frame means patients can switch cushion types without the need to be resized for an entirely new mask. This can not only help reduce out-of-pocket costs for patients but can help DMEs run their practice with less inventory and carrying costs.



Patients found Philips Therapy Mask System 3100 NC and 3100 SP less intimidating than their prescribed mask^{1,2}



¹Analysis after 10 days of use of July 2020 Nasal Patient Preference trial data n=123. Prescribed masks include Wisp, Activa, AirFit N10, AirFit N20, ComfortGel Blue Nasal, EasyLife, Eson, Eson 2, Mirage FX, Pico, TrueBlue.

²Analysis after 10 days of use of July 2020 Pillows Patient Preference trial data n=113. Prescribed masks include AirFit P10, Brevida, Nuance/Nuance Pro, Opus, Pilairo, and Swift/Swift FX/Swift LT.

³Data analysis of June 2020 Clinician Ease of Use and First Impression trial where n=31. Comparison masks were the ResMed AirFit P10 and the ResMed AirFit N30.

⁴Data analysis of June 2020 Clinician Ease of Use and First Impression trial – optional blindfold task where n=27. Comparison masks were the ResMed AirFit P10 and the ResMed AirFit N30

⁵Based on engineering analysis of the size and weight of the Philips Therapy Mask System 3100 NC and 3100 SP, AirFit N30 and AirFit P10 masks.

⁶Per the Design History File, Philips Therapy Mask 3100 NC and 3100 SP were designed with specific features and benefit goals including correlated sizing between cushion types, comfort, minimal size, adjustable headgear, sound, and ease of use.

⁷When compared to leading CPAP mask brands

Philips Therapy Mask System 3100 NC and 3100 SP

Lightweight design
Our smallest, lightest
adult mask ever¹

Open field of vision⁵
Wear glasses, read
and watch TV

Silicone Pillows
Small, stable, and quiet⁴

**Only tube-in-front
mask system²**
Easily switch
between cushions

Integrated magnets
Allow patients to easily
attach cushions^{3,4}


¹Vs. Philips adult nasal and pillows masks. Based on engineering analysis of size of Philips' adult nasal and pillows masks available in the US as of April 1, 2020.

²When compared to leading CPAP mask brands

³Analysis after 10 days of use of July 2020 Nasal Patient Preference trial data n=123. Prescribed masks include Wisp, Activa, Airfit N10, Airfit N20, ComfortGel Blue Nasal, EasyLife, Eson, Eson 2, Mirage FX, Pico, TrueBlue.

⁴Analysis after 10 days of use of July 2020 Pillows Patient Preference trial data n = 113. Prescribed masks include Airfit P10, Brevida, Nuance/Nuance Pro, Opus, Pilairo, and Swift/Swift FX/Swift LT.

⁵Design of the mask prevents the cushion from touching the bridge of the nose

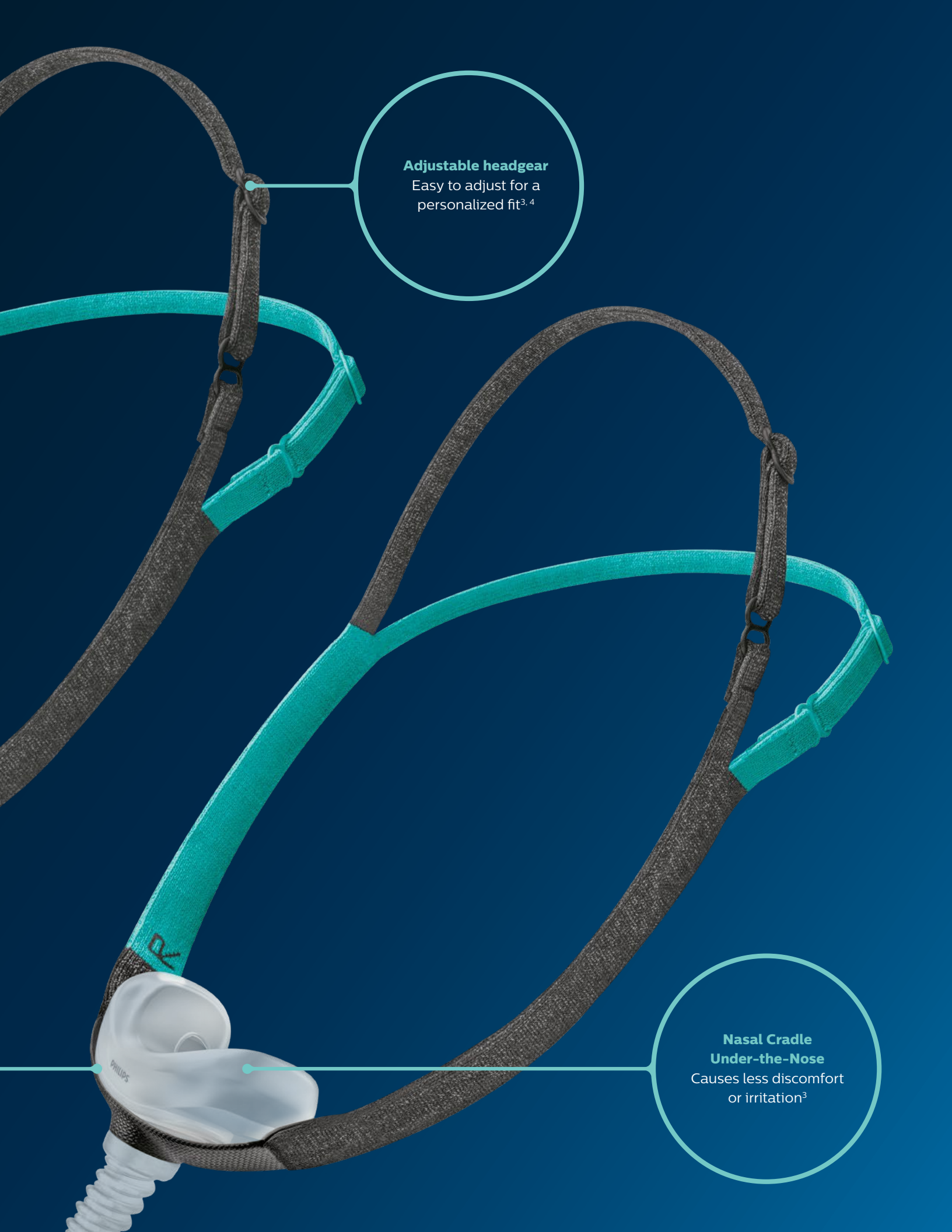

Adjustable headgear
Easy to adjust for a personalized fit^{3,4}

Nasal Cradle Under-the-Nose
Causes less discomfort or irritation³

Two cushion types, one mask frame

Philips Therapy Mask System 3100 NC and 3100 SP is an innovative tube-in-front CPAP mask system that lets you and your patients switch between a nasal cradle under-the-nose or silicone pillows cushion while keeping the same frame. Innovative, polarized magnets in the cushion and frame allow effortless assembly of the cushions to the frame.¹ When patients feel the magnets click into place, they can be confident their cushion is correctly attached. The intuitive design with orientation markers on the headgear for placement on the head can minimize patient mask errors that lead to callbacks, saving you time.



The Silicone Pillows cushion, developed using data from 3D facial-scanning technology, are designed to accommodate various nostril shapes.



The Nasal Cradle Under-the-Nose cushion features a minimal design that eliminates soreness across the nose bridge.²





Clinicians say Philips Therapy Mask 3100 SP is easier to use than the leading pillows brand³

¹Data analysis of June 2020 Clinician Ease of Use and First Impression trial – optional blindfold task where n=27. Comparison masks were the ResMed AirFit P10 and the ResMed AirFit N30

²Design of the mask prevents the cushion from touching the bridge of the nose.
³Among n=31 RTs and STs, vs ResMed AirFit N30 and P10. Data analysis of June 2020 Clinician Ease of Use and First Impression trial where n=31. Comparison of masks where the ResMed AirFit P10 and the ResMed AirFit N30.



Philips Therapy Mask System 3100 NC and 3100 SP patient-fitting guide

Follow the instructions to get started fitting Philips Therapy Mask System 3100 NC and 3100 SP

Nasal Cradle Under-the-Nose

Silicone Pillows

 <p>1</p>	<p>For the best performance, use the smallest cushion that fits the nose. Hold the cushion in front of the face. Check that the L and R symbols are upright and on the correct sides.</p>	
 <p>2</p>	<p>Holding the cushion against the nose with one hand, pull the headgear straps over the head.</p>	
 <p>3</p>	<p>The cushion opening will be directly under the nostrils. At no time should the nose be inside the cushion opening.</p>	
 <p>4</p>	<p>Use the sliders to adjust the top headgear straps. Pull the slider to either the left or right to tighten or loosen.</p>	
 <p>5</p>	<p>Adjust the bottom headgear straps.</p>	
 <p>6</p>	<p>The mask should fit loosely and comfortably against the face.</p>	

Ordering information

Global Philips Therapy Mask System 3100 NC and 3100 SP fit-pack

	Part number		
Philips Therapy Mask 3100 NC, fit-pack, w/HGR, GBL	1144608		
Philips Therapy Mask 3100 SP, fit-pack, w/HGR, GBL	1144610		
Philips Therapy Mask 3100 NC with headgear	Part number	Philips Therapy Mask 3100 SP with headgear	Part number
Extra small, Philips Therapy Mask 3100 NC, w/HGR, GBL	1145459	Extra small, Philips Therapy Mask 3100 SP, w/HGR, GBL	1145464
Small, Philips Therapy Mask 3100 NC, w/HGR, GBL	1145460	Small, Philips Therapy Mask 3100 SP, w/HGR, GBL	1145465
Medium, Philips Therapy Mask 3100 NC, w/HGR, GBL	1145461	Medium, Philips Therapy Mask 3100 SP, w/HGR, GBL	1145466
Large, Philips Therapy Mask 3100 NC, w/HGR, GBL	1145463	Large, Philips Therapy Mask 3100 SP, w/HGR, GBL	1145478
Medium wide, Philips Therapy Mask 3100 NC, w/HGR, GBL	1145462	Medium wide, Philips Therapy Mask 3100 SP, w/HGR, GBL	1145467
Philips Therapy Mask 3100 NC without headgear	Part number	Philips Therapy Mask 3100 SP without headgear	Part number
Extra small, Philips Therapy Mask 3100 NC, w/oHGR, GBL	1145044	Extra small, Philips Therapy Mask 3100 SP, w/oHGR, GBL	1145055
Small, Philips Therapy Mask 3100 NC, w/oHGR, GBL	1145045	Small, Philips Therapy Mask 3100 SP, w/oHGR, GBL	1145056
Medium, Philips Therapy Mask 3100 NC, w/oHGR, GBL	1145046	Medium, Philips Therapy Mask 3100 SP, w/oHGR, GBL	1145057
Large, Philips Therapy Mask 3100 NC, w/oHGR, GBL	1145047	Large, Philips Therapy Mask 3100 SP, w/oHGR, GBL	1145058
Medium wide, Philips Therapy Mask 3100 NC, w/oHGR, GBL	1145048	Medium wide, Philips Therapy Mask 3100 SP, w/oHGR, GBL	1145059
Philips Therapy Mask 3100 NC replacement cushion	Part number	Philips Therapy Mask 3100 SP replacement cushion	Part number
Extra small, Philips Therapy Mask 3100 NC, cushion, RP	1145064	Extra small, Philips Therapy Mask 3100 SP, cushion, RP	1145069
Small, Philips Therapy Mask 3100 NC, cushion, RP	1145065	Small, Philips Therapy Mask 3100 SP, cushion, RP	1145070
Medium, Philips Therapy Mask 3100 NC, cushion, RP	1145066	Medium, Philips Therapy Mask 3100 SP, cushion, RP	1145071
Large, Philips Therapy Mask 3100 NC, cushion, RP	1145067	Large, Philips Therapy Mask 3100 SP, cushion, RP	1145072
Medium wide, Philips Therapy Mask 3100 NC, cushion, RP	1145068	Medium wide, Philips Therapy Mask 3100 SP, cushion, RP	1145073
Philips Therapy Mask 3100 NC replacement parts	Part number	Philips Therapy Mask 3100 SP replacement parts	Part number
Philips Therapy Mask System 3100 NC-SP, frame, RP	1145074	Philips Therapy Mask System 3100 NC-SP, frame, RP	1145074
Philips Therapy Mask System 3100 NC-SP, headgear, RP	1145075	Philips Therapy Mask System 3100 NC-SP, headgear, RP	1145075
Philips Therapy Mask System 3100 NC-SP, swivel, RP	1145076	Philips Therapy Mask System 3100 NC-SP, swivel, RP	1145076
Philips Therapy Mask System 3100 NC-SP, sizing gauge, RP	1145077	Philips Therapy Mask System 3100 NC-SP, sizing gauge, RP	1145077
International Philips Therapy Mask 3100 NC fit-pack	Part number	International Philips Therapy Mask 3100 SP fit-pack	Part number
Philips Therapy Mask 3100 NC, fit-pack, w/HGR, INT	1145038	Philips Therapy Mask 3100 SP, fit-pack, w/HGR, INT	1145049
International Philips Therapy Mask 3100 NC replacement parts	Part number	International Philips Therapy Mask 3100 SP replacement parts	Part number
Philips Therapy Mask 3100 NC, fit-pack, w/HGR, CE	1145039	Philips Therapy Mask 3100 SP, fit-pack, w/HGR, CE	1145050
Philips Therapy Mask 3100 NC, fit-pack, w/HGR, RC	1145040	Philips Therapy Mask 3100 SP, fit-pack, w/HGR, RC	1145051
Philips Therapy Mask 3100 NC, fit-pack, w/HGR, AP	1145041	Philips Therapy Mask 3100 SP, fit-pack, w/HGR, AP	1145052
China, Japan and South Korea Philips Therapy Mask 3100 NC fit-pack	Part number	China, Japan and South Korea Philips Therapy Mask 3100 SP fit-pack	Part number
Therapy Mask 3100 NC, fitpack, w/HGR, CN	1145042	Therapy Mask 3100 SP, fitpack, w/HGR, CN	1145053
Therapy Mask 3100 NC, fitpack, w/HGR, JP	1145043	Therapy Mask 3100 SP, fitpack, w/HGR, JP	1145054
Therapy Mask 3100 NC, fitpack, w/HGR, KR	1149794	Therapy Mask 3100 SP, fitpack, w/HGR, KR	1149795



Learn more

For more information, visit
philips.com/professional/3100-NC-SP.



© 2021 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.

www.philips.com/respironics

Caution: U.S. federal law restricts these devices to
sale by or on the order of a physician.

RRDPGH MK 10/14/21 PN 1148609 MCI 4110223 v02
1010 Murry Ridge Lane, Murrysville, PA 15668
800 345 6443 • 724 387 4000